Deta
Elis
Holding

Gesundheit
Innovationen
Mensch

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DeVita Combi

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Languages

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Complete Recovery
All-In-One Device

DeVita Combi

User Manual • Warranty Document

Purpose of the Device

The all-in-one device for complete recovery is designed to eliminate the strain caused by various intruders in the human body: viruses, bacteria, fungi, protozoa, and helminths. It also shields against stressors, maintains homeostasis, enhances performance, and combats premature aging.

The device works in 5 different fields of action

DeVita AP

The operation of the AP device is rooted in the pioneering research of engineer R. Rife and subsequent advancements by his successors. Each living organism possesses a unique electromagnetic signature. The device reproduces waves at specific frequencies, corresponding to the inherent frequency of target microorganisms. These waves penetrate the body, resonate with the microorganisms' frequencies to amplify their natural oscillations. The resulting powerful oscillations have detrimental effects on pathogenic microorganism and parasites. The external AP unit emits waves at precise frequencies tailored to various types of viruses, bacteria, helminths, protozoa, and fungi, facilitating thorough detoxification of the body.

DeVita RITM

The operation of the Ritm device based upon extensive research conducted by R. Voll, H. Schimmel, AF Kramer, and other distinguished scientists. Their findings indicate that each healthy organ has a distinct frequency, which differs from that of a weakened counterpart. Consequently, the device emits low-intensity frequencies tailored to various organs and systems, mirroring those of healthy organs. These frequencies effectively "reprogram" the body towards a harmonic range of healthy frequencies, thereby enhancing adaptive reserves and restoring self-regulatory processes.

QUANT

The device's impact on the body is based on research within the field of quantum theory of light. Visible light comprises a wide spectrum of colors, each with unique wavelengths and frequencies. These colors can affect specific organs or systems through skin sensors, based on the principles of light interaction with biological tissues.

HV electrotherapy

The device operates by utilizing electrical pulses or impulses, offering HV STEP and HV RAMPING electrotherapy modes for the treatment of various pathological conditions. These operating options are integrated directly into the device for ease of use.

Caution! The device may only be used by children under adult supervision. Children are not allowed to play with the device. Persons with physical or mental disabilities may use the device only under the supervision of appropriate personnel.

Package contents


All-in-one device "DeVita COMBI"..... 1 pc.
External AP unit (anti-parasitic action fig. 1)..... 1 pc.
External RITM unit (regulatory action fig. 2)..... 1 pc.
External QUANT unit (fig. 3)..... 1 pc.
Hand electrodes (fig. 5)..... 2 pcs.
Foot electrodes (fig. 6)..... 2 pcs.
Cable set (fig. 7)..... 1 pc.
External Power supply unit (fig. 4)..... 1 pc.
Bag for carrying and storing the device..... 1 pc.
User Manual / Warranty Document (passport)..... 1 pc.


Fig. 1 Ext. AP unit

Fig. 2 Ext. Ritm unit

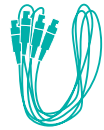1

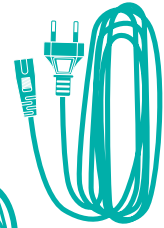2

Fig. 3 Ext. QUANT unit

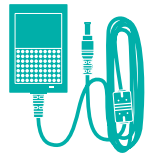
Fig. 4 Power supply unit

3

4

Fig. 7 Cable set

Fig. 5 Hand electrodes

Fig. 6 Foot electrodes

Technical Specifications

External AP unit

Output waveform Sinusoidal

Frequency range, Hz 0,1 – 72000

Maximum frequency adjustment error, not more than, % ±0.01

Carrier frequency, MHz 27

Overall dimensions, not more than, mm
Diameter 110

Weight, not more than, Kg 0,1

External RITM unit

Output waveform Sinusoidal

Frequency range, Hz 0,1 – 1000

Maximum frequency adjustment error, not more than, % ±0.01

Carrier frequency, kHz 10

Overall dimensions, not more than, mm
Diameter 110

Weight, not more than, Kg 0,1

QUANT unit

Variation of repetitive pulses, Hz 0,1 – 1000

Maximum frequency adjustment error, not more than, % ±0.01

Optical range, MHz 468,525,590, 660, 880

Power supply:

- battery Li-Pol 3.7V, pcs 1
- continuous battery operation time, not more than, hours 10
- maximum battery charge time, hours 8

Overall dimensions, not more than, mm 95 x 45 x 25

Weight, not more than, Kg 0,1

HV STEP and HV RAMPING Electrotherapy

Shape of the output signal: pulses square positive or negative and pulses square positive or negative with hyperbolic attenuation

Maximum pulse width, not less than, V.....200
Pulse duration at the level of 0.5, ms.....1
Frequency range, Hz.....0,1 - 100
Maximum frequency adjustment error, not more than, %±0.01

Power supply unit:

AC - DC adapter.....GSM18B05-P1J 220V\50Hz, 5V\3A

Device “DeVita Combi”

Keyboard Type.....touch sensor
Screen Type.....TFT 240x320
Overall dimensions, not more than, mm..... 240x185x104
Weight:
without peripherals and accessories, not more than, kg.....0,8 in
Complete set, not more than, kg.....2,5

Operating conditions:

Ambient temperature.....from 10 to 35˚C
Relative humidity.....80% at a temperature of 25˚C
Shelf life, not less than, years.....5

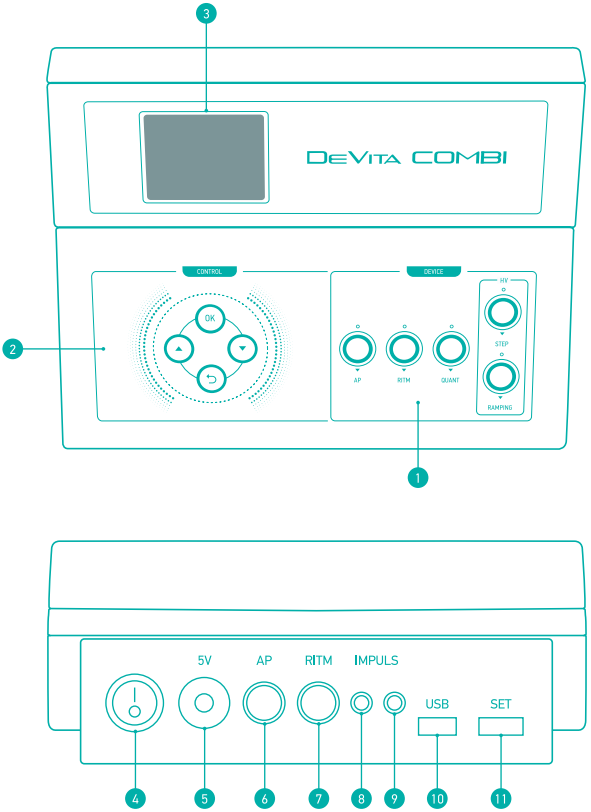
Fig. 8 External appearance of the device:
front view (top), back view (bottom)

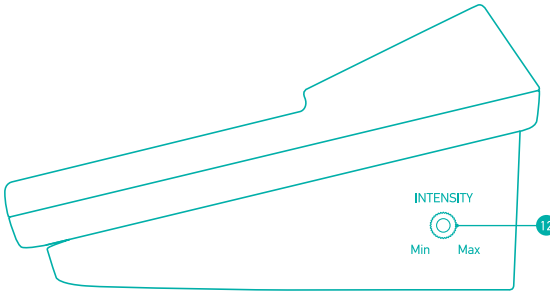
Fig. 9 External appearance of the device
(right side view)

Device composition

1. “DEVICE” Action Selection Panel
2. “CONTROL” panel
3. Display
4. Power ON/OFF switch
5. External power supply connection port (“DC 5V “)
6. External AP module unit connection port
7. External RITM unit connection port
8. HV electrode connection port “-” (black)
9. HV electrode connection port “+” (red)
10. Micro-USB port for connection to a computer.
11. USB type A — for setting-up and charging the QUANT unit
12. INTENSITY — control knob to adjust the pulse amplitude of HV unit.

Terms of use

1. Do not wipe the surface of the device with alcohol and other organic solvents, remove the dirt with light movements with a soft brush or a cotton swab slightly dampened with water.

2. Avoid positioning the device in close proximity to excessively hot, dusty, or humid environments.

3. Following storage in a cold room or transition to winter conditions, allow the device to acclimate to room temperature for at least one hour before attempting to activate it.

Preparing the device for operation

1. Connect the power supply to the "5V", power port (No. 5 in Fig. 8) of the device. With cable No. 2, connect the power supply to a 220V outlet.


2. Press the ON-OFF switch (No. 4 in Fig. 8), while the indicator LED of one of the modes (actions) and the display (No. 3 in Fig. 8), will light up.


3. Connect the AP and/or RITM external units to the device with cables No.1. Connect the electrodes (hand or foot) to cables No. 4. Cables numbers are shown in Figure 7.

4. Set the INTENSITY control knob to the MIN position.

Using the device


After activating the device, the screen will display information about the type of action selected and the corresponding function indicator will light up on the “DEVICE” panel. To change the action unit, use the action selector buttons on the “DEVICE” panel.

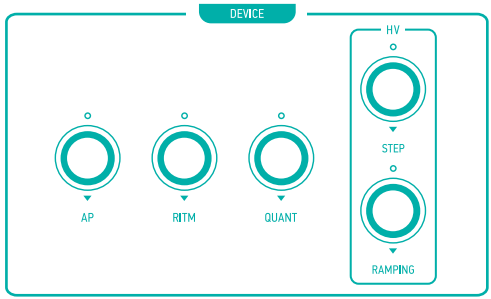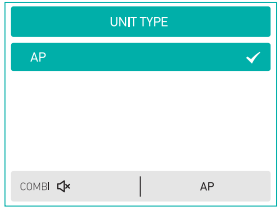
The bottom of the screen includes information about the volume status as well as information about the type of action unit selected.

Working with the external AP unit

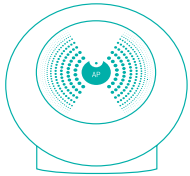
The external “AP” unit is connected and controlled via a cable. Connect the external “AP” unit with cable No. 1 to the “AP” port on the back of the “DeVita Combi” device (see Fig. 8).

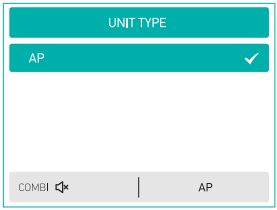Press the ‘AP’ button on the “DEVICE” panel.

OK

Press the ‘OK’ button on the ‘CONTROL’ panel to continue.

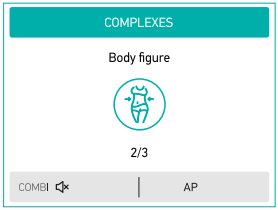
The SET of programs selection window will appear.

The window displays the current program SET name, its icon, and its number. The total number of remaining SET programs is also displayed.

To select a SET of programs, use the

“CONTROL” panel navigation buttons.

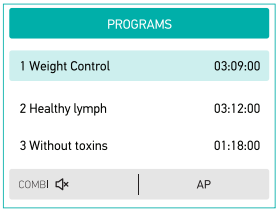
When you press the “OK” button, the program selection window will appear on the screen.

Program selection is made with the navigation buttons on the “CONTROL” panel.

When the “OK” button is pressed, the program starts.

Before starting, you need to make sure that the external AP is connected.

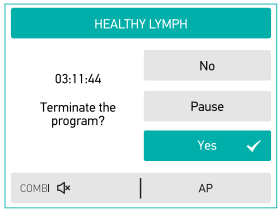The program launch window contains the name of the program, the time remaining until its completion and a menu for managing the program's operation.

When the program starts, the LED on the AP external unit should light up indicating that the unit is connected correctly.

By selecting "Pause" and pressing the "OK" button, the program is temporarily stopped. In this case, the display shows the time until the end of the pause.

The pause is stopped by selecting the "Yes" option and pressing the button "OK" and then the program will continue to run.

If pause is enabled and "Yes" is not selected, then the program stops when the pause time elapses. The device switches to the program selection window.


Working with the external RITM unit

The external "RITM" unit is connected and controlled via a cable.


Connect the external "RITM" unit with Cable No. 1 in the "RITM" port on the back of the "DeVita Combi" device (see Fig. 8).

Press the “RITM” button on the “DEVICE” panel.

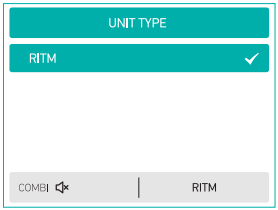
Working with the external “RITM” unit is similar to working with the external “AP” unit.


Working with the QUANT unit

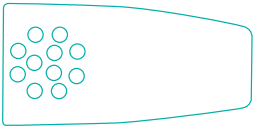Operation with the external QUANT unit is carried out via a Bluetooth wireless channel. Ensure that the QUANT unit number matches the device number. To start, press the “QUANT” button on the “DEVICE” panel.

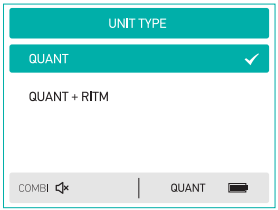Operation with the “QUANT” unit can occur in two modes: independently or in conjunction with the “RITM” unit. Navigate using the buttons on the “CONTROL” panel to select your desired mode. The battery status of the external “QUANT” unit is indicated in the lower right corner for easy monitoring.

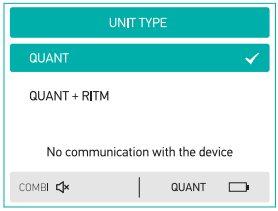If there is no connection via Bluetooth, the message “No connection to the device” appears on the screen and the battery image changes.

Without a connection, accessing the SET program selection becomes impossible.

To choose SET programs or select programs, follow the same procedure as described for the external “AP” unit.

When the “QUANT+RITM” mode is selected, ensure the connection of the external “RITM” unit is verified.

Working with the HV unit


The device includes an electrical pulse generator, located inside the ‘DeVita Combi’ device, and external hand and foot electrodes.


The electrodes are connected with the cables No. 4 to the HV (+red, -black) connectors on the back side of the device (see Fig. 8).

The “INTENSITY” control knob, located on the lower right side of the device (see Fig. 9), adjusts the amplitude of the pulses.

The “HV” unit has the capability to operate independently or in conjunction with the external “RITM” unit. In such instances, the “HV” unit offers two operational modes: “HV STEP” (with stepping pulses) and “HV RAMPING” (with accelerating pulses).

To switch between the “HV STEP” and “HV RAMPING” modes, simply press the corresponding button on the “DEVICE” panel.

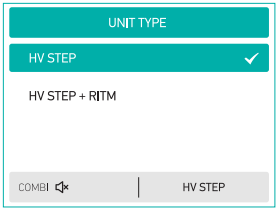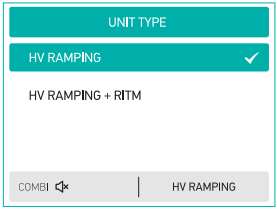
The selection of SET programs and programs is carried out similarly to the method described for the external "AP" unit.

Settings Section


To go to the settings menu, press the "back" button.

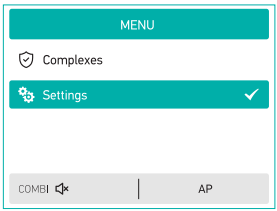Use the navigation buttons to select "Settings" and press "OK" on the "CONTROL" panel.


The settings selection menu will appear on the screen.

Language selection

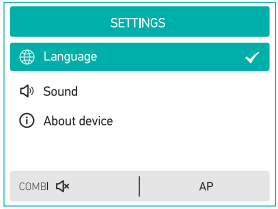
Use the navigation buttons to select ‘Language’. Click the ‘OK’ button.

The language settings menu appears on the screen.

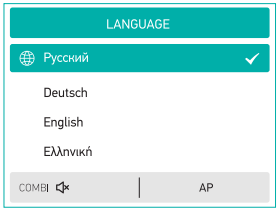
Use the navigation buttons to select the desired language. Click the ‘OK’ button.

Sound Setting

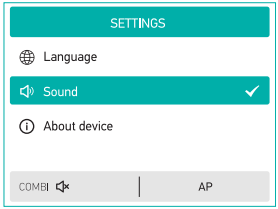
Use the navigation buttons in the settings menu to select “Sound”. Press the “Okay” button.

The “Sound” menu appears on the screen.

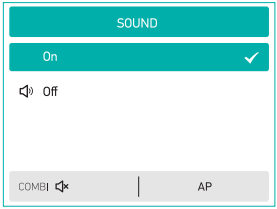Use the navigation buttons to select “On” or “Off”. Click the ‘OK’ button. Selecting “Off” the sound of the keys will be muted.


About the device

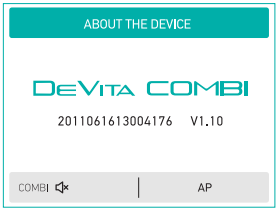Use the navigation buttons to select “About the Device”. Click the ‘OK’ button.

Your device build-in number and software version will appear on the screen.


Settings menu when using the QUANT unit

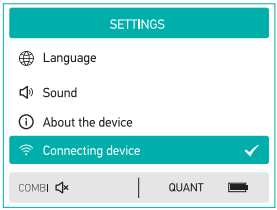
When working with the QUANT unit, an additional “Connect” option appears in the “Settings” menu.

The option is intended to connect a new QUANT unit in case of loss or damage of the existing one.

To install, connect the QUANT unit to the “SET” port on the back side of the “DeVita Combi” with cable No. 3

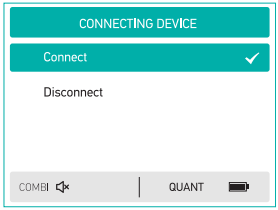
On the device screen, use the navigation buttons to select the "Connect" option Click the 'OK' button.

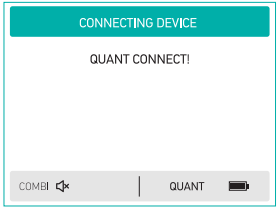If the connection is successful, the message "QUANT CONNECTED!" should appear on the screen.

(Connection of more than 1 unit is not allowed).

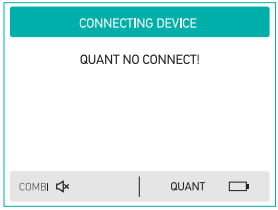To disconnect an unused QUANT unit, follow the procedure outlined above by connecting it to the device and accessing the "Disconnect" option.


Upon completion of the disconnection process, the message "QUANT NOT CONNECTED!" will be displayed on the screen, confirming the successful disconnection.

Charging the QUANT unit

The QUANT unit is equipped with a battery, and its charge level can be assessed through the charge indicators located at the bottom of the display.

It is important to note that operation with the device is not possible when the charging indicator is red.

Charging the battery is facilitated by connecting the device to a computer or an external charger using cable No. 3.

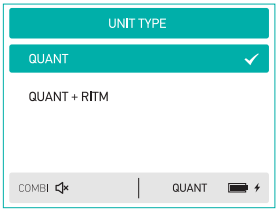At the same time, a charging icon in the form of a "lightning bolt" will appear to the right of the charging indicator.

Safety information

Device is:

- non-toxic;

- not flammable;

- non-explosive;

When charging the battery, Class 3 devices must be powered exclusively by an external low-voltage AC-DC adapter.

Environmental Protection

The built-in battery of the QUANT unit contains substances that may pollute the environment. At the end of the device's life cycle and before handing it in for recycling, the battery must be removed. The battery must be disposed of at designated collection points. 

Safety instructions

CAUTION: Never connect the device to a power supply (adapter) if the power cord is damaged. Avoid attempting to repair the device independently. In the event of malfunction, please contact your local technical service centers for assistance.

CAUTION: Direct exposure of the QUANT unit to the eyes may cause discomfort.

CAUTION: To avoid accidental high voltage exposure of the electrotherapy when switching on the current, it is necessary to set the "INTENSITY" control knob to the MIN position.

Save

The device, when removed from the package, must be stored in a room with an air temperature of 10 to 35 °C and a relative humidity of no more than 80%.

To avoid interruptions, it is recommended not to leave the power supply unit connected to the device. It is recommended that you keep the packaging and all manuals of the device during the warranty period.

Transporting the device

Since the device contains a display that is sensitive to external mechanical influences, it is recommended during transport:

- to protect the device from drops
- to protect the device from sudden shocks
- to make sure that no other objects fall on the device

The device must be protected from moisture and the effects of active chemicals.

Manufacturer's Warranty

The manufacturer guarantees the compliance of the "DeVita Combi" all-in-one wellness device, with the technical specifications, provided that the consumer follow the rules of operation, transport, and storage.

The warranty period of the device is 24 months from the date of purchase from a retail store.

If the date of sale and the seller's seal are absent in the manual, the warranty period is calculated from the date of production by the manufacturer.

During the warranty period, owners are entitled to free repairs upon presentation of the warranty repair form. Warranty repairs are conducted exclusively at the manufacturer's facilities. However, the cost of shipping the defective device is the responsibility of the buyer.

Please note that repairs covered by warranty will not be accepted without the submission of an application for warranty repair and a certificate of acceptance. Additionally, violations of the device's seals may void warranty coverage.

The warranty does not cover damages resulting from:

- Poor handling;
- Natural wear and tear;
- Breach of security seals;
- External defects such as cracks or bumps.

The buyer has the right to replace the defective device with a new one in the following cases:

- the device has been repaired three times during the warranty period;
- the device cannot be repaired and warranty repair forms are included.

For technical information contact exclusive distributor:

Data Elis Cyprus

Address: 179 Agias Fylaxeos St, 1st floor, office 12, Limassol, 3882 Tel: +357 25222910

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EC Declaration of Conformity
Nr.W17

CE

The company **DETA-ELIS EUROPA GMBH** with a seat at Kaiserstr. 74, 63065 Offenbach, Germany

Herewith to DECLARE

That the equipment described below complies with all applicable at the time of construction essential requirements of the **Directive 2011/65/EU (RoHS)** and the **Directive 2014/30/EU (EMC).**
This declaration relates exclusively to the equipment in the state in which it was placed on the market, and excludes alterations or manipulations, components which are added and/or operations carried out subsequently by the final user without our official consent.

Equipment: DeVita
Type / Model: Combi
Standards Applied: IEC 61000-4-2-2008
GOST P МЭК 60601-2:2014
GOST 52770:2007
GOST ISO 10993-1-2-5-10-11-12-18-2011

Date: 01 October 2018

Authorized Representative: The General Manager Ms. Konopleva Tatyana

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Deta-Elis Europa GmbH
Deta-Elis Europa GmbH Kaiserstrasse 74,